

Contents


General Tips	2
How to make FaSSIF	3
How to make FeSSIF	4
How to make FaSSGF	5
How to make FaSSIF-V2	6
How to make FeSSIF-V2	7
How to make Dog FaSSIF	8
How to make Dog FaSSGF	9
How to make FaSSCoF	10
How to make FeSSCoF	11

General Tips



Do not refrigerate or freeze the media

It's important to use clean glassware for making the media to prevent microbial growth

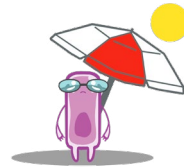

Keep the media out of direct sunlight and preferably in the dark

Store Buffer Concentrate bottles at room temperature and ensure both cap and nozzle are tightly closed after use



Once opened, store the bottle of powder in the refrigerator (2-8 °C) with the cap tightly closed



FaSSIF



STATE
Fasted

FLUID SIMULATED
Small Intestinal

CONTAINS
Bile Salts, Lecithin

BUFFER
Phosphate

DATE OF BIRTH
1998

COMPOSITION

Component	Concentration (mM)
Taurocholate	3
Phospholipids	0.75
Sodium	148
Chloride	106
Phosphate	29



How to make 1 L of FaSSIF

For other quantities of FaSSIF and instructions how to prepare your own buffer, scan this QR code or go to biorelevant.com/FFF



Prepare buffer



- Weigh 41.65 g (or 35.00 mL) of FaSSIF Buffer Concentrate into suitable container
- Make up to volume (1.000 L) with purified water at room temperature

Check pH

- If required adjust the pH to 6.5 with Sodium hydroxide 1M or Hydrochloric acid 1M




Add powder



- Weigh 2.240 g of FaSSIF/FeSSIF/FaSSGF powder into suitable container
- Add buffer to volume (1.000L)
- Stir until powder is completely dissolved

Ready to use

- Let stand for 2 hours. It will become slightly opalescent. Your FaSSIF is now ready to use
- We recommend using the medium at room temperature or at 37°C. It can be used for up to 48 hours at these temperatures

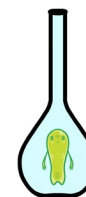

FeSSIF



STATE
Fed

FLUID SIMULATED
Small Intestinal

CONTAINS
Bile Salts, Lecithin

BUFFER
Acetate

DATE OF BIRTH
1998

COMPOSITION

Component	Concentration (mM)
Taurocholate	15
Phospholipids	3.75
Sodium	319
Chloride	203
Acetic acid	144



How to make 1 L of FeSSIF

For other quantities of FeSSIF and instructions how to prepare your own buffer, scan this QR code or go to biorelevant.com/FFF



Prepare buffer



- Weigh 81.41 g (or 70.00 mL) of FeSSIF Buffer Concentrate into suitable container
- Make up to volume (1.000 L) with purified water at room temperature

Check pH

- If required adjust the pH to 5.0 with Sodium hydroxide 1M or Hydrochloric acid 1M



Add powder



- Weigh 11.20 g of FaSSIF/FeSSIF/FaSSGF powder into suitable container
- Add buffer to volume (1.000L)
- Stir until powder is completely dissolved

Ready to use

- FeSSIF will have a clear appearance and is now ready to use
- We recommend using the medium at room temperature or at 37°C. It can be used for up to 48 hours at these temperatures



FaSSGF



STATE
Fasted

FLUID SIMULATED
Stomach

CONTAINS
Bile Salts, Lecithin

BUFFER
None (uses Hydrochloric Acid)

DATE OF BIRTH
2005

COMPOSITION

Component	Concentration (mM)
Taurocholate	0.08
Phospholipids	0.02
Sodium	34
Chloride	59



How to make 1 L of FaSSGF

For other quantities of FaSSGF and instructions how to prepare your own buffer, scan this QR code or go to biorelevant.com/FFF



Prepare buffer



- Weigh 36.78 g (or 35.00 mL) of FaSSGF Buffer Concentrate into suitable container
- Make up to volume (1.000 L) with purified water at room temperature

Check pH

- If required adjust the pH to 1.6 with Hydrochloric acid 1M



Add powder



- Weigh 0.0597 g of FaSSIF/FeSSIF/FaSSGF powder into suitable container
- Add buffer to volume (1.000L)
- Stir until powder is completely dissolved

Ready to use

- FaSSGF will have a clear appearance and is ready to use
- We recommend using the medium at room temperature or at 37°C. It can be used for up to 48 hours at these temperatures



FaSSIF-V2



STATE
Fasted

FLUID SIMULATED
Small Intestinal

CONTAINS
Bile Salts, Lecithin

BUFFER
Maleate

DATE OF BIRTH
2008

COMPOSITION

Component	Concentration (mM)
Taurocholate	3
Phospholipids	0.2
Sodium	106
Chloride	69
Maleic acid	19



How to make 1 L of FaSSIF-V2

To calculate other quantities of FaSSIF-V2 scan this QR code or visit biorelevant.com/V2FAS



Prepare buffer

In about 0.9 L of purified water dissolve:

- 1.392 g of Sodium hydroxide pellets (NaOH)
- 2.220 g of Maleic acid (C₄H₄O₄)
- 4.010 g of Sodium chloride (NaCl)



Check pH

- Adjust the pH to 6.5 with Sodium hydroxide 1M or Hydrochloric acid 1M
- Make up to volume (1 L) with purified water at room temperature




Add powder

- Weigh 1.790 g of FaSSIF-V2 powder into suitable container
- Add buffer to volume (1.000L)
- Stir until powder is completely dissolved



Ready to use

- Let stand for 1 hour. Your FaSSIF-V2 is now ready to use
- We recommend using the medium at room temperature or at 37°C. It can be used for up to 48 hours at these temperatures

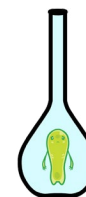

FeSSIF-V2



STATE
Fed

FLUID SIMULATED
Small Intestinal

CONTAINS
Bile Salts, Lecithin, Digestive Components

BUFFER
Maleate

DATE OF BIRTH
2008

COMPOSITION

Component	Concentration (mM)
Taurocholate	10
Phospholipids	2
Oleate	0.8
Glycerol monoleate	5
Sodium	218
Chloride	125
Maleic Acid	55



How to make 1 L of FeSSIF-V2

For other quantities of FeSSIF-V2, scan this QR code or go to biorelevant.com/V2FES



Prepare buffer

In about 0.9 L of purified water dissolve:

- 3.270 g of Sodium hydroxide pellets (NaOH)
- 6.390 g of Maleic acid (C₄H₄O₄)
- 7.330 g of Sodium chloride (NaCl)



Check pH

- Adjust the pH to 5.8 with Sodium hydroxide 1M or Hydrochloric acid 1M
- Make up to volume (1 L) with purified water at room temperature



Add powder

- Weigh 9.760 g of FeSSIF-V2 powder into suitable container
- Add buffer to volume (1.000L)
- Stir until powder is completely dissolved



Ready to use

- Let stand for 1 hour. FeSSIF-V2 will have a clear appearance and is now ready to use
- We recommend using the medium at room temperature or 37 °C. It can be used for up to 48 hours at these temperatures



Dog FaSSIF



STATE
Fasted

FLUID SIMULATED
Small Intestinal (Canine)

CONTAINS
**Bile Salts, Lecithin, Digestive
Components**

BUFFER
Phosphate

DATE OF BIRTH
2013

COMPOSITION

Component	Concentration (mM)
Taurocholate	5
Taurodeoxycholate	5
Phospholipids	1.25
Lysophospholipids	1.25
Oleate	1.25
Sodium	121
Chloride	60
Phosphate	29



How to make 1 L of Dog FaSSIF

To calculate other quantities of Dog FaSSIF scan this QR code or visit biorelevant.com/DOGFAS



Prepare buffer

In about 0.9 L of purified water dissolve:

- 0.870 g of Sodium hydroxide pellets (NaOH)
- 3.950 g of Monobasic sodium phosphate monohydrate ($\text{NaH}_2\text{PO}_4 \cdot \text{H}_2\text{O}$)
- 3.480 g of Sodium chloride (NaCl)



Check pH

- Adjust the pH to 7.5 with Sodium hydroxide 1M or Hydrochloric acid 1M
- Make up to volume (1 L) with purified water at room temperature




Add powder

- Weigh 7.800 g of Dog FaSSIF/Dog FaSSGF powder into suitable container
- Add buffer to volume (1.000L)
- Stir until powder is completely dissolved



Ready to use

- Let stand for 1 hour. Your Dog FaSSIF is now ready to use
- We recommend using the medium at room temperature or at 37°C. It can be used for up to 48 hours at these temperatures

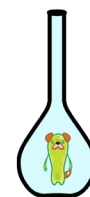

Dog FaSSGF



STATE
Fasted

FLUID SIMULATED
Stomach (Canine)

CONTAINS
Bile Salts, Lecithin, Digestive
Components

BUFFER
Maleate

DATE OF BIRTH
2013

COMPOSITION

Component	Concentration (mM)
Taurocholate	0.1
Taurodeoxycholate	0.1
Phospholipids	0.025
Lysophospholipids	0.025
Oleate	0.025
Sodium	59
Chloride	15
Maleic	22



How to make 1 L of Dog FaSSGF

For other quantities of Dog FaSSGF and instructions how to prepare the pH 1.5 version, scan this QR code or go to biorelevant.com/DOGFAS



Prepare buffer

In about 0.9 L of purified water dissolve:

- 1.610 g of Sodium hydroxide pellets (NaOH)
- 2.520 g of Maleic acid (C₄H₄O₄)
- 1.100 g of Sodium chloride (NaCl)



Check pH

- Adjust the pH to 6.5 with Sodium hydroxide 1M or Hydrochloric acid 1M
- Make up to volume (1 L) with purified water at room temperature



Add powder

- Weigh 0.160 g of Dog FaSSIF/Dog FaSSGF powder into suitable container
- Add buffer to volume (1.000L)
- Stir until powder is completely dissolved




Ready to use

- Your Dog FaSSGF will have a clear appearance and is ready to use
- We recommend using the medium at room temperature or at 37°C. It can be used for up to 48 hours at these temperatures



FaSSCoF

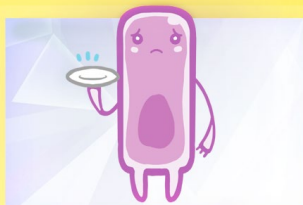

STATE
Fasted

FLUID SIMULATED
Colonic

CONTAINS
**Bile Salts, Lecithin, Digestive
Components**

BUFFER
Tris/Maleate

DATE OF BIRTH
2010

COMPOSITION

Component	Concentration (mM)
Sodium Cholate	0.15
Phospholipids	0.3
Oleate	0.1
Sodium hydroxide	120
TRIS	45
Maleate	76



How to make 1 L of FaSSCoF

To calculate other quantities of FaSSCoF visit this QR code or visit biorelevant.com/COFAS



Prepare buffer

In about 0.9 L of purified water dissolve:

- 4.800 g of Sodium hydroxide pellets (NaOH)
- 8.800 g of Maleic acid ($C_4H_4O_4$)
- 5.500 g of Tris base ($C_4H_{11}NO_3$)



Check pH

- Adjust the pH to 7.8 with Sodium hydroxide 1M or Hydrochloric acid 1M
- Make up to volume (1 L) with purified water at room temperature




Add powder

- Weigh 0.340 g of FaSSCoF powder into suitable container
- Add buffer to volume (1.000L)
- Stir until powder is completely dissolved



Ready to use

- Let stand for 1 hour. Your FaSSCoF is now ready to use
- We recommend using the medium at room temperature or at 37°C. It can be used for up to 24 hours at these temperatures

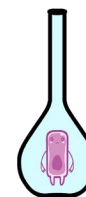

The FaSSCoF powder is very hygroscopic so it's important to replace the bottle cap immediately after use and keep it tightly closed



FeSSCoF



STATE
Fed

FLUID SIMULATED
Colonic

CONTAINS
**Bile Salts, Lecithin, Digestive
Components**

BUFFER
Tris/Maleate

DATE OF BIRTH
2010

COMPOSITION

Component	Concentration (mM)
Sodium Cholate	0.6
Phospholipids	0.5
Oleate	0.2
Sodium hydroxide	34
TRIS	31
Maleate	30



How to make 1 L of FeSSCoF

To calculate other quantities of FeSSCoF scan this QR code or visit biorelevant.com/COFES



Prepare buffer

In about 0.9 L of purified water dissolve:

- 0.660 g of Sodium hydroxide pellets (NaOH)
- 3.500 g of Maleic acid ($C_4H_4O_4$)
- 3.690 g of Tris base ($C_4H_{11}NO_3$)



Check pH

- Adjust the pH to 6.0 with Sodium hydroxide 1M or Hydrochloric acid 1M
- Make up to volume (1 L) with purified water at room temperature



Add powder

- Weigh 0.740 g of FeSSCoF powder into suitable container
- Add buffer to volume (1.000L)
- Stir until powder is completely dissolved



Ready to use

- Let stand for 1 hour. Your FeSSCoF is now ready to use
- We recommend using the medium at room temperature or at 37°C. It can be used for up to 24 hours at these temperatures



The FeSSCoF powder is very hygroscopic so it's important to replace the bottle cap immediately after use and keep it tightly closed





biorelevant.com

Suite 4

The QMB Innovation Centre

42 New Road

London

E1 2AX

United Kingdom

Tel: +44 (0)20 7790 5328

Email: info@biorelevant.com

